premium

D&D UM-200/UM-300 premium

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● Product Description

Cordless Rechargeable Massager 01

Cordless, rechargeable and portable massager. It is designed for massage of any body parts according to chosen massage sticks

Convenient Long-lasting Lithium-ion batteries 02

Powerful lithium-ion batteries generate the same level of massage strength compared to electricity-powered massagers

Durable, Outstanding Motor 03

Smaller, lighter motor, but powerful performance

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● More information about the Battery

- Batteries equipped : 7.2V Lithium-ion 2200mA 18650 2 Cell
- Battery life varies by ways of application and/or operating conditions.
- If the product is not in use for a long time, please keep its power off and check the battery on a regular basis to avoid possible natural discharge. Batteries last longer this way.
- Batteries don't stand fire or humidity. Please keep them away from heat and avoid using in places where the surrounding temperature is relatively hot.
- Keep the power off and store the product in a cool place.

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● Application / Recharging methods

1. Depending on where to massage on the body, attach the massage stick of one's choice.
2. Green light is on as the power on.
Set the speed and start massaging. Please be aware that the recharger needs to be pulled off from the product, otherwise it won't work.
3. For User's safety, the product is designed to automatically power off and stops itself if too much loads are applied on the product. In this case, power off first - if it is still on - plug in the adapter to an electrical outlet and then plug in the recharger on the product. Wait for about three seconds and then pull off both, and start massaging.

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
● Application / Recharging methods

- Battery life : continuous use for 180 minutes (if no load is applied)
- Battery tests were done in the manufacturer's facility with its own equipment.
- Battery life varies depending on ways of application and/or the surrounding conditions

4. Plug in the adapter provided by the manufacturer to an electricity outlet, pull off the charging terminal cap, and then insert the charger. The red light indicates the charging is in process.
5. The green light indicates the recharging is complete. Remove the charger and put the cap back into the charger input.

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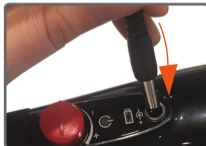
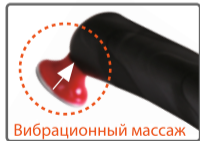
● Vibration Massage Use (Applicable only with CM-07)

1. Turn the dial and set it at '0' speed level. The red light is on immediately and the beep sound will hear after 3 seconds.
2. The beep sound indicates the vibration module at the end of the hand-grip is activated. Set the dial at your preferred speed.
3. Skin Massage Stick, Scalp Massage Stick, and Facial Massage stick are available.

-
- Excessive load on the product during massage should activate the safety module and shut down the product. In this case, power off and back on. If the product sti doesn't work, plug in the charger and plug it off after 2 or 3 seconds.
 - If the product is kept in a cool place and not in use for a while, please don't massage too heavy for about 10 minutes, for a sudden motor stop can take place.
-

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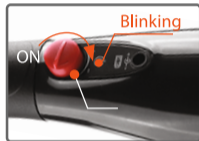● Safety warning

For your safety and to avoid any property damage. please read this manual carefully and use as indicated.

- Please use the recharger provided by the manufacturer. Unspecified adapter or unlicensed power generators might cause explosion, electric shock, or fire.
- Please avoid using the product while driving, for it might cause a traffic accident.
- Please avoid using the massager with wet hands or inserting metallic material into the charging socket for it might cause an electric shock.
- Please keep the product away from heating appliances such as oven, microwave, hot cooking utensils, or high pressure container, for it might cause a leak in the batteries and / or product malfunction.
- Please do not apply heavy forces, take apart, or drill the product or batteries, for it can cause an electric shock, explosion or fire.
- Please stop using the product and contact our customer service immediately for the following cases: unusual smell, or noise, or smoke from the product or batteries.
- Please avoid using the product on a heat-generating mat, for it can overheat the batteries and cause malfunction or explosion.
- Please keep the product out of reach of children or pets.
- The recommended temperature to use the product are between 5°C (41°F) and 35°C (95°F). Please keep the product in places where the temperature is not too cold nor too hot, otherwise it can cause malfunction or explosion.

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● *Safety warning*

- Do not drop, apply excessive forces or put something heavy over the product, for it can cause malfunction and/or damages.
- Do not take apart, modify or change the parts for it can cause an electric shock and the product is no longer subject to Warranty.
- Do not use the product other than the usage stated in the manual, for it can cause breakdown.
- Do not keep or use the product in where temperature changes drastically or where humidity is high, for it can cause damage in electrical circuits and batteries.
- Do not keep the product near high magnetic field or where it can be affected by magnetic field, for it can cause inconvenience from discharged batteries.
- If the product is not in use for a while, check the batteries once every month and recharge them. If the batteries remain discharged too long, recharge cycle and performance might decrease and it can cause breakdown, explosion and/or fire.
- Use the product appropriately in public.
- If the product is in use commercially for 20 minutes or over, it might be a cause for a breakdown.
- It is recommended not to use the product during pregnancy.

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● Warranty

Our warranty covers only those defects that arise as a result of normal use of the product and does not cover any other problems, including those that arise as a result of following cases:

1. Improper or inadequate maintenance or modification
2. Accident, misuse, abuse, contamination or other external causes
3. The use of unspecified adapter and accessories
4. Loss or damage in transit
5. Damages that occurs as a result of one's failure to follow the instructions

This Warranty does not apply to expendable or consumable parts and does not extend to any product from which the serial number has been removed.

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Name of product		Warranty terms
Model name	D&D UM-200/UM-300 premium	One year
Manufacturing No.		
Date purchased		
Customer signature		

Warning

AC Charger specification

Input voltage : 100V~240V 50/60Hz 0.4A Max

Output voltage : 8.5V DC 1A

Time to recharge : 90 minutes

Warning : Please use the adapter provided by the manufacturer. Time to recharge the batteries varies by the surrounding conditions and usually it takes about 90 minutes.

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