*premium*

D&D UM-100 premium

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Cordless Rechargeable Massager

01

Cordless Rechargeable Massager: Fast recharging portable percussion massager. It is designed for massage of any body parts except head and face

Convenient Long-lasting Lithium-ion battery

02

Small and light but powerful: Efficiently designed percussion motor to be lightweight yet yields outstanding performance and power.

Durable, Outstanding Motor

03

Self-stand design and three massage tips: Without the need of another composition of stand and with three various massage tips.

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
## ● Battery Specification



- Battery: 3.7V Lithium-ion 2600mA 18650 1 Cell.
- Battery life: continuous use up to 120 minutes (if no load is applied).
- Battery tests were done in the manufacturer's facility with its own equipment. Battery life varies depending on usage.
- If the product is not in use for a long time, please keep its power off. Batteries will last longer if kept recharged versus at a discharged state.
- Please keep the massager away from heat and avoid using in places where the surrounding temperature is relatively hot.
- During storage keep the power off and store the product in a cool place.

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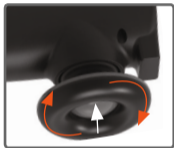
## ● Application

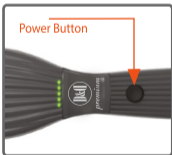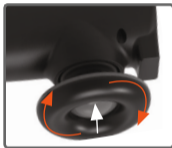1. Depending on where to massage on the body, attach the massage stick of one's choice.
2. Green light is on as the power on (Please know that there are about two seconds of waiting time) Set the speed and start massaging. Please be aware that the charger needs to be pulled off from the product, otherwise it won't work.
3. For User's safety, the product is designed to automatically power off and stops itself if too much loads are applied on the product. In this case, power off first- if it is still on - plug in the charger to an electrical outlet and then plug in the charger on the product. Wait for about three seconds and then pull off charger, and start massaging.
4. Plug in the adapter provided by the manufacturer to an electricity outlet, and then insert the charger. The red light on the charger indicates the charging is in process.
5. The green light on the charger indicates the recharging is complete. Remove the charger.

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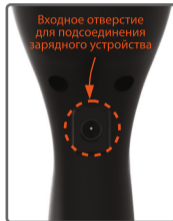## ● Charging and Using Body Elements

For the best results charge massager before use.

1. Plug the CHARGER into an electrical outlet.
2. Plug the CHARGER into the CM-03 Massager.
3. The charger has an LED light which will display red, orange or green to indicate charge state of massager. A green light on the charger indicates 100% charge.
4. Remove CHARGER for operation.
5. Push the On/Off SPEED BUTTON for a second until green light is on. Set desired speed by pushing On/Off SPEED BUTTON which varies by 4 steps.
6. The 3 attachments screw on clockwise to the unit. They should be screwed on hand tight. Do not forcefully over tighten.
7. For user's safety, the product may power off if too much load is applied. In this case, Pure-Wave may need to be reset by plugging it into the charger for 3 seconds.

- 
- Battery life : continuous use for 120 minutes (if no load is applied)
  - Battery life varies depending on ways of use and/or the surrounding conditions.
- 

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## ● *Warnings and Safety*

**For your safety and to avoid any property damage. please read this manual carefully and use as indicated.**

- Only use the charger provided. Unspecified charger or adapter may cause overheating, explosion, electric shock or fire.
- Only use certified Body elements massage stick adaptors for your safety and to prevent damage of Body elements cordless massager.
- Do not use the product while driving a vehicle .
- Do not use the massager with wet hands. Do not insert any object into the charger input. This may casue damage to the product or electric shock.
- Keep the product away from heating appliances such as oven, microwave, hot cooking utensils, high pressure cooker. This may cause a degradation of the battery and/or product malfunction.
- Do not apply heavy force to the product.
- Do not disassemble or drill into the product and/or battery. This may cause product malfunction, electric shock or fire.
- Do not use product on heat generating material. This may cause the battery to overheat and cause malfunction or overheating.
- Keep product out of reach of children or pets.
- Do not drop, apply excessive force or heavy weight over Body elements massager for it can cause malfunction and/or damages.
- The recommended temperature to use the product in places where the temperature is not too cold nor too hot. 5°C (41°F) to 35°C (95°F). Use common sense when using or storing the Body elements. Leaving in overly hot places (such as inside a parked car on hot day or next to a window with direct sun) may cause malfunction or overheating to battery.

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## ● *Safety warning*

- Do not use massager on:
  - a. Sensitive skin, swollen or inflamed areas, in the presence of poor circulation, where skin eruptions are present, or in the presence of unexplained calf or abdominal pain.
  - b. Baricose veins
  - c. The throat area
  - d. Genital areas
  - e. Unconscious or sleeping person
- Do not take apart, modify or change the parts. This may cause electric shock and the product will no longer be covered under manufacturer warranty.
- Do not use the product other than the usage stated in the manual, for it can cause breakdown.
- We recommend not to use this as a sexual personal massager. This is a powerful massager and may cause severe bodily injuries.
- Do not keep or use the product in areas where temperature changes drastically or where humidity is high. This may cause damage to electrical circuits and battery.
- If the product is not in use for extended duration, charge the unit once every month. If the battery remain discharged too long, this may degrade the life cycle and performance of product.
- Do not keep product near a high magnetic field or where it can be affected by a magnetic field, for it can cause battery to discharge.
- This product is intended for consumer home use. Commercial use will cause premature breakdown of the product.

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## ● Warranty (1-Year Limited)

---

**Our warranty covers only those defects that arise as a result of normal use of the product and dose not cover any other problems, including those that arise as a result of following cases :**

---

1. Improper or inadequate maintenance or modification
2. Accident, misuse, abuse, contamination or other external causes
3. The use of unspecified adapter and accessories
4. Loss or damage in transit
5. Damages that occurs as a result of one's failure to follow the instructions

This Warranty does not apply to expandable or consumable parts and does not extend to any product from which the serial number has been removed.

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Name of product		Warranty terms
Model name	D&D UM-100 premium	One year
Manufacturing No.		
Date purchased		
Customer signature		

## Warning

### AC Charger specification

**Input voltage** : 100V~240V 50/60Hz 0.3A

**Output voltage** : 5V DC 1A

**Time to recharge** : 120 minutes

**Warning** : Please use the adapter provided by the manufacturer. Time to recharge the batteries varies by the surrounding conditions and usually it takes about 120 minutes.

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